


Operating & Maintenance Manual

Gas Control Panel

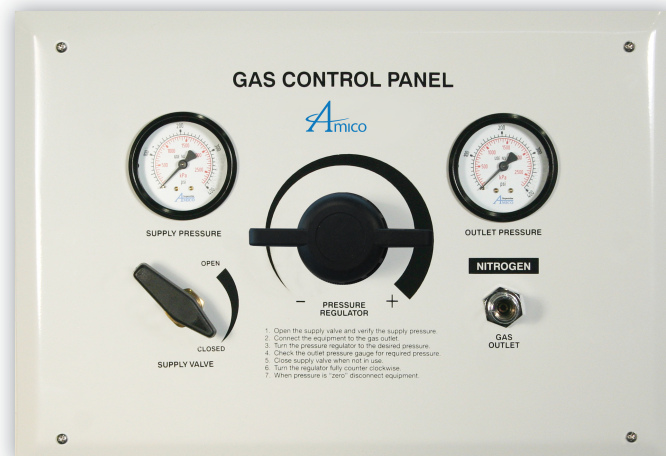

Contents

Intended Use/Introduction	3
Operation	4
Installation	4
Parts List	4 - 5
Notes	6 - 7

Intended Use

The Amico Gas Control Panel is designed to deliver Nitrogen or Air to turbo-surgical tools.

Introduction

The Panel consists of a Shut-off valve, Delivery and Supply pressure Gauges, Pressure Regulator and DISS Outlet, in addition an outlet pipe is supplied that can be piped to additional Gas outlets.

The panel comes in two styles "Relieving" and "Non-relieving". A relieving panel allows gas to vent from the regulator, to reduce the pressure from the OUTLET PRESSURE side of the panel. CSA requires that a non-relieving panel be used where CSA Standards are practiced.

Operation

1. Turn the Shut-off valve to the “Open” position.
2. Verify that the correct supply pressure is available on the supply gauge.
3. Attach a tool to the DISS outlet and tighten the nut (do not over tighten) or HP “Safety-Swing Coupling” outlet.
4. Adjust the pressure regulator to the required pressure by observing the Outlet Gauge.
5. Operate the surgical tool and observe the outlet pressure and adjust if needed.
6. When the control panel is not in use reduce the outlet pressure to “zero” and turn off the supply valve.


Installation

Mount the back-box into the wall and adjust the adjustable side mounting brackets so that the top of the box will be flush with the finished drywall. Connect the inlet pipe to the gas pipeline and the outlet to additional outlets (this is optional). The pipeline can now be pressure tested to the required testing pressure. The outlet High Pressure Hose is shipped with one end connected to the outlet and the other end connected to the inlet on the nitrogen back box. This will enable the installer to pressure test the complete pipeline, including remote outlets, without having to attach the front assembly.

To attach the control panel front assembly to the control panel back box you must first disconnect the braided hose from the male PT adapter that is attached to the inlet pipe on the control panel back box. Then loosen and remove the panel cap from the male PT adapter that is marked “output” on the back of the control panel front assembly. Attach the braided hose from the back box to the male PT adapter (avoid over-tightening). Connect the “input” braided hose from the back of the control panel front assembly to the “input” of the male PT adapter in the back box (avoid over-tightening). Mount the control panel front assembly to the back box with the enclosed machine screws.

Parts List

Nitrogen control panel English DISS (Relieving type)	N-CONP-E-REL
Nitrogen control panel English DISS (Non-Relieving type)	N-CONP-E
Nitrogen control panel English “Dyna-Quip” HP (Relieving type)	N-CONP-E-HP-REL
Nitrogen control panel English “Dyna-Quip” HP (Non-Relieving type)	N-CONP-E-HP
Nitrogen control panel French DISS (Non-Relieving type)	N-CONP-F
Instrument Air control panel English DISS (Non-Relieving type)	N-CONP-E-I
Instrument Air control panel English DISS (Relieving type)	N-CONP-E-REL-I

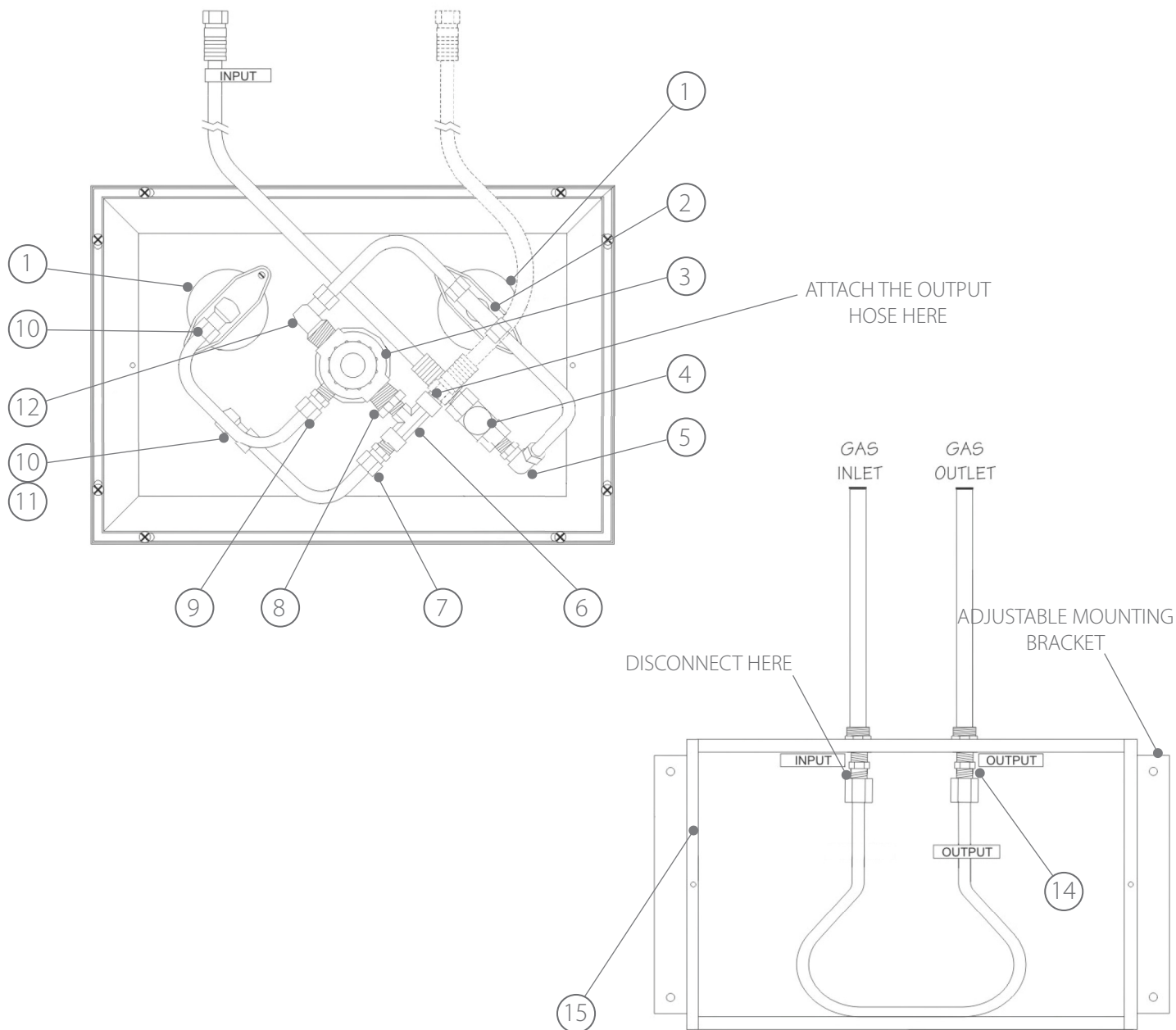

Item	Description	Part Number	Item	Description	Part Number
1	Gauge 0-400 psi	M-X-MAN-52B	9	Connector 3/8"	N-X-ELB-ALM2508
2	Tee (female)	N-X-TEE-AL-F304	10	Elbow 1/4"	N-X-ELB-ALF2508
3a.	Regulator Relieving	N-X-REG-W-HP-R	11a.	DISS Check - Nitrogen	S-DIS-DEMC-NIT
3b.	Regulator Non-Relieving	N-X-REG-W-HP	11b.	DISS Check - Instrument Air	S-DIS-DEMC-IAIR
4	Shut-off Valve	N-X-VALVE	11c.	Dyna-Quip Check	N-X-CON-ADA-B25
5	Brass Elbow	N-X-ELBW-AL-M37	12	Elbow 1/2" * 3/8"	N-X-ELBW-AL-305
6	Tee (female)	N-X-TEE	13	Braided Pigtail	N-X-PIGTAIL
7	Connector 3/8"	N-X-CON-AL-375	14	Pigtail/Pipe Adapter	N-X-PT-ADAPTER
8	Brushing 1/2" * 3/8"	N-X-BUSH	15	Gas Panel Back Box	N-X-BOX

Notes

Notes

www.amico.com

Amico Pipeline | www.amico.com

85 Fulton Way, Richmond Hill
Ontario, L4B 2N4, Canada

Toll Free Tel: 1.877.264.2697

Toll Free Fax: 1.866.440.4986

Tel: 905.764.0800

Fax: 905.764.0862

